

Mobile Endoscopic Inspection Table

Belflit Compact




Concept



Belflit seamlessly supports the entire endoscopy process, including sedation. Everything from pre-op preparation and inspection/treatment to patient transport to the recovery room—all with just this single table.

Layout Sample

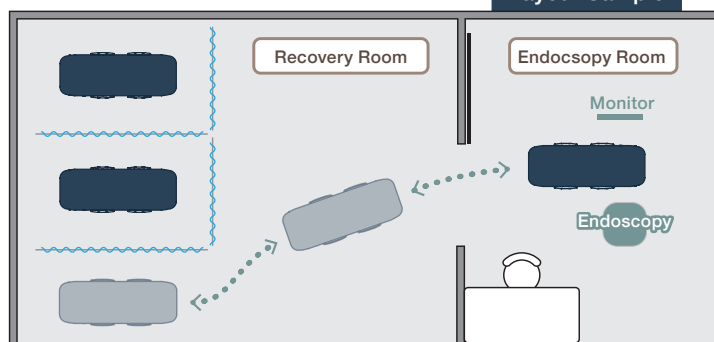

Compact Footprint Design

Our compact design saves space in standby position, ensures smooth patient entry and wheelchair transfers, and delivers effortless maneuverability.



Quick Reclining Operation

Effortless operation supported by a gas cylinder regardless of a patient's physical build. The slip-resistant mechanism that lowers the seat section as the backrest is raised prevents the patient from sliding.



Wide Height Range

A 430 mm height range (500–930 mm) ensures low height for easy patient access and the optimal working height for clinicians.



Ample Space for Legroom

The compact table base maximizes floor space, allowing closer patient access and hassle-free placement of multiple foot controllers.




Comfortable Mattress

Provides soft and comfortable cushioning

OPTION



A profiled mattress made of specially designed urethane foam with an uneven surface can be selected.

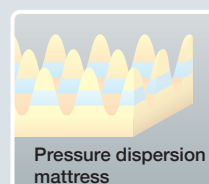

Quick and Effortless Setup

It is designed to support emergency transport and side-by-side stretcher positioning.


One-Handed Operation

Each fence can be operated with one hand. When moved to the use or storage position, it is automatically locked by releasing it. (It cannot be unlocked even when pushed from inside.)



Rotatable Side Fence

The rotatable side fences do not require space when switching between the use and storage positions. When not in use, they can be stored below the tabletop, which provides easy accessibility during treatments or when transferring the patient onto or off of the table.

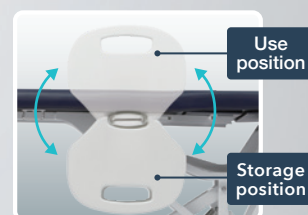

Better Mobility with Durable Casters

Single-pedal operation enables switching between four-wheel free, directional-travel, and four-wheel lock.

When totally locked, casters do not rotate or swivel, ensuring stability during examination and treatments.



Lithium-ion Battery

It uses a compact lithium-ion battery with a long runtime and durability for repeated recharging compared to lead-acid batteries.

The table can be moved up and down approximately 100 times on a single charge (when new).

The audible alarm and LED indicator will alert you when the battery level is low.



Hand Controller with Two Memory Settings

With two preset switches, you can adjust the table to your preferred height with a single press.



OPTION



Foot Controller with Two Memory Settings

The table can also be operated with an optional foot controller.

Choice of Layouts to Suit Your Operations

■ Caster Lock Pedal


In addition to the standard equipment (two), additional installations are available as an option. Mounting positions can be selected from the four locations shown in the figure on the right, A - D.



OPTION

■ IV POLE

Up to 2 locations can be selected from the installation locations 1 to 4 in the figure on the right. In addition, the position 5 or 6 can be selected to follow the elevation of the main unit.

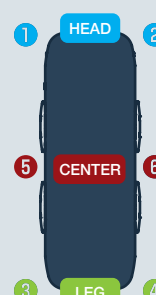


*Combination of 1 and 2 or 3 and 4 is not available.


Options



SIDE FENCE COVER
EB-BFT-SC

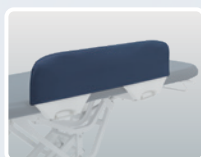

HEADREST CUSHION
EB-BFT-M

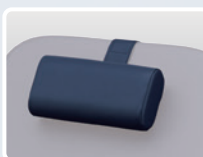

FOOT CONTROLLER
EB-BFT-FS



CASTER LOCK PEDAL



LITHIUM-ION BATTERY
DR-LIBP





BATTERY CHARGER FOR LITHIUM-ION BATTERY
DR-LICHG



PRESSURE DISPERSION MATTRESS

A profiled mattress with pressure dispersion effect can be selected.

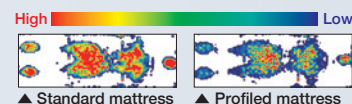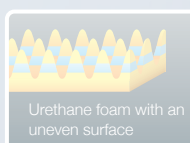

Pressure redistribution test result (In-house test)
Testee : Male (165cm, 67kg)

TRAY
EB-BFTC-KM



OXYGEN CYLINDER HOLDER
EB-BFTC-BH

* Holder Dimensions
-Length : 440 mm
-Inside Diameter : 109 mm
* Oxygen cylinder is not included.



IV POLE (Bracket 1pc is included.)

HEAD

EB-BFT-DN

BRACKET FOR IV POLE
EB-BFT-DN-BK



CENTER

EB-BFTC-DC

BRACKET FOR IV POLE
EB-BFTC-DC-BK



LEG

EB-BFTC-DL

BRACKET FOR IV POLE
EB-BFTC-DL-BK



Color Selection



ES23
ALDO BLUE



ES11
OLIVE BROWN



ES15
LIGHT GREEN



ES16
PURPLE BLUE

*Choose from the four color options shown above.

Antibacterial performance and Alcohol resistance



Antibacterial performance

Tested in accordance with JIS Z 2801:2012
Antimicrobial products—Test for antimicrobial activity and efficacy (Test bacteria: Staphylococcus aureus and Escherichia coli)



Alcohol resistance

Tested in accordance with our standard protocol

Specifications

Model	EB-BFTC
Dimensions	650mm×1850mm (tabletop)
Table height range	500mm-930mm (Electric motor)
Side fence	Rotatable & retractable
Caster	Four-wheel free, Directional-travel, Total lock
Electrical rating	DC21.6V 130VA
Product weight	85kg
Standard accessories	Hand controller, Battery pack, Charger

Dimensions

