

# Doctor Stool

Functional and durable doctor stools professionally support surgery and consultation. Simple and easy foot operation for height adjustment keeps hands clean.

## Hydraulic Pump Stool




Height adjustment can be controlled from the foot pedal while seated.

### LS-150HYKB



#### Seat and backrest angle adjustment

The backrest and seat can be tilted to the desired angle by operating the lever. Tilting the seat raises the position of the lower back, which takes pressure off the back of the thighs. The bucket-shaped seat helps maintain a stable posture.

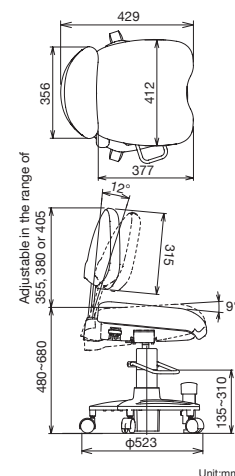

OPTION



#### Armrest LS-OPT-AR

Adjustable height and forward/backward position. Cushions gently support the user's arms.

### LS-150HY



OPTION



#### Backrest LS-OPT-SB

Rotating bar-shaped backrest.




#### Round Backrest LS-OPT-MB

Rounded backrest for secure back support (Does not rotate.)



#### Armrest LS-OPT-AR

Adjustable height and forward/backward position. Cushions gently support the user's arms.

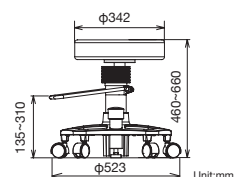

## Gas Cylinder Stool




The chair can be lowered by depressing the ring at the base and raised by taking your weight off the chair while depressing the ring.

### LS-009GS



OPTION


#### Backrest LS-OPT-SB

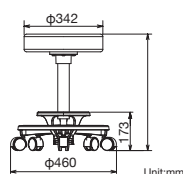

### LS-007GSB



A height-adjustable backrest is included as a standard accessory.

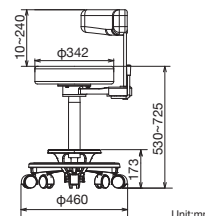

Model	LS-150HYKB	LS-150HY	LS-009GS	LS-007GSB
Base	Hydraulic pump	Hydraulic pump	Gas cylinder	Gas cylinder
Seat height	480mm-680mm	460mm-660mm	505mm-700mm	530mm-725mm
Weight	22kg	15kg	8kg	13kg
Backrest	Standard	Optional LS-OPT-SB / LS-OPT-MB	Optional LS-OPT-SB	Standard
Armrest	Optional LS-OPT-AR	Optional LS-OPT-AR	-	-