


Examination Table

# Flalith Series



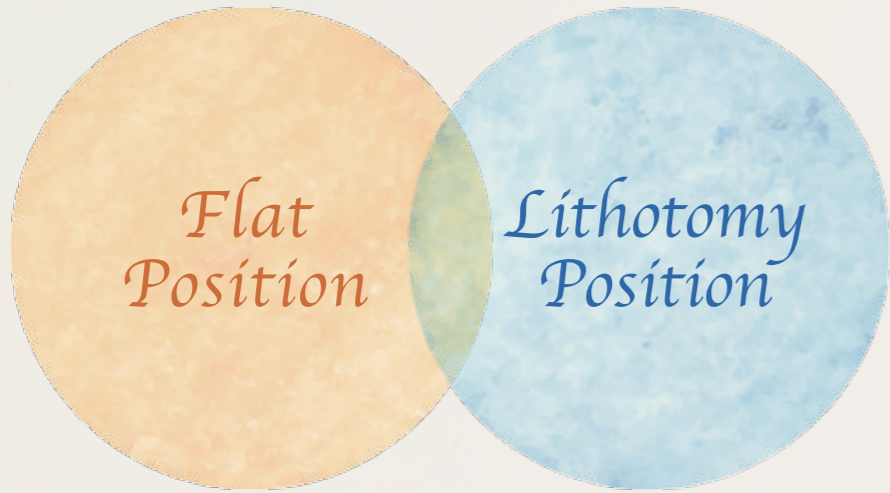

In your Obstetrics & Gynecological care, if the clinical examination and abdominal ultrasound are performed in the consultation room, and internal examination and transvaginal ultrasound are performed in the internal examination room, we recommend the type of outpatient care that is completed in one room, enabling both clinical and internal examinations to be performed on one examination table.




# Flalith


“Flalith” is a word coined by TAKARA BELMONT combining “flat position” and “lithotomy position.” We aimed at ease of use in these two positions.



Layout examples for reference

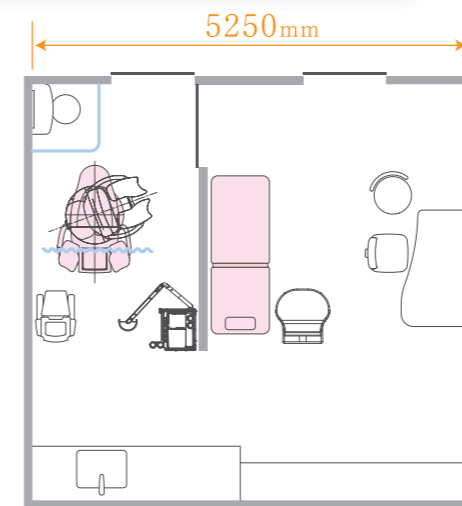

[Clinical and internal examinations separated]

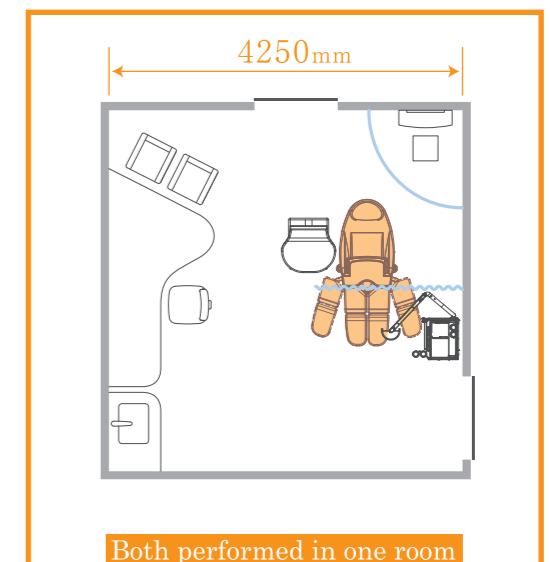

Both performed in one room

# Flalith DG-FL750

Being able to perform both clinical and internal examinations on one table, you can use the space and time effectively. The soft form of the table also supports your practice by helping relax your patients.

**Seating position**  
Low initial height makes it easy to get on and off for patients.



**Lithotomy Position**  
Switched with one touch to the posture in clinical examination/transvaginal ultrasound inspection.



### Wider examination area

Easy to use transvaginal probe



### Selective leg opening

For patients whose leg/hip joint has limitation to move outward, each leg support can open individually and gradually from the totally closed state.



### Higher pelvis position

The seat part can be tilted independently. This can help make the treatment field more secure for doctor.



**Flat Position**  
Switched with one touch to the posture in clinical examination/abdominal ultrasound inspection.



Simple side design makes it easy to get closer to your patient

End part of the leg supports bends to make it easy to raise knees when opening legs, enabling natural opening of legs.

## Function

### Wireless foot control

With the wireless foot control using infrared rays, the area around the base of the examination table is stress-free.



### Main unit switch

- Abdominal ultrasound position
- Internal examination position
- Standby position
- Selective leg opening



### Plastic tray

Corrosion-resistant and water-repellent. Detachable allowing for maintenance.

# Flalith *K* DG-FL750K

The lifting-type leg supports that can facilitate various examinations and procedures.

Urological treatment and procedures can also be performed in the table.

## Lithotomy Position

Securing a wide space for examination allowing easy use of cystoscope.



Can be switched from the seating position to each examination position by one touch of the foot control.

## Flat Position

Low initial height makes it easy to get close to your patients and facilitates examinations.



### Function

#### Wireless foot control

With the wireless foot control using infrared rays, the area around the base of the examination table is stress-free.



#### Higher pelvis position

The seat part can be tilted independently. This can help make the treatment field more secure for doctor.

#### Selective leg opening

For patients whose leg/hip joint has limitation to move outward, each leg support can open individually and gradually from the totally closed state.


### Option

#### Arm board



DG-FL750-A

#### Side fence

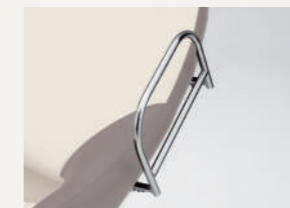

DG-FL750-SF

#### Disposable roll sheets



DG-OPT-FLR

#### Drain bucket



DG-OPT-DB (5L/10L)

### Color Variation



ES9 [ Ivory ]






ES23 [ Aldo blue ]






ES20 [ Coral pink ]

### Specifications

Flalith  
DG-FL750N

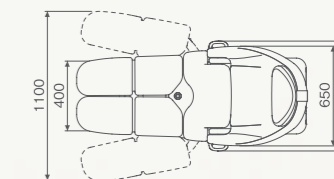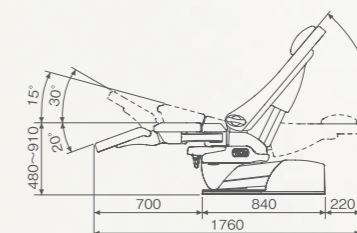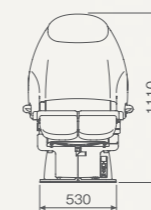

Flalith *K*  
DG-FL750K

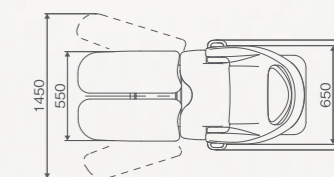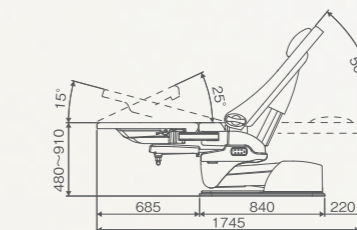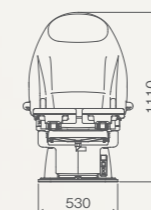

( mm )

Table height range	480mm – 910mm
Back support angle	0°- 50°
Seat angle	0°- 15°
Preset/Auto-return operation	○
Rating	100 VAC/110 V/120 V, 50/60 Hz, 5 A    220 VAC /230 V/240 V, 50/60 Hz, 3 A
Net weight	Flalith:205kg / Flalith-K:208kg

## Disposable Examination Skirt

Use of a special skirt for examination reduces the burden on patients from pulling up and down their clothes. The wraparound skirt is easy for pregnant patients to wear.

Length : 65cm  
Size : Free (abdominal circumference ~120cm)  
Material : Polypropylene





<https://www.takarabelmont.co.jp/english/medical/index.html>

**TAKARA BELMONT CORPORATION**

2-1-1, Higashishinsaibashi, Chuo-ku, Osaka, 542-0083, Japan  
TEL +81-6-6213-5945 FAX +81-6-6212-3680

Book no. M-017  
Printed in Japan 2020. 4