

DELIVERY TABLE DG-880

Highly functional delivery table equipped with enhanced functions capable of supporting emergency deliveries in its relaxing and healing shape



Leather Color Selection



■ P1. Pink



■ B1. Brown

EMERGENCY DELIVERIES

Highly efficient functions capable of accommodating emergencies

60cm Wide Tabletop

Easy to approach for operation

- Table top is 60cm wide with the side fences removed.
- The table height is adjustable from 58cm to 102cm.
- Arm boards can be easily attached/detached in a single operation.

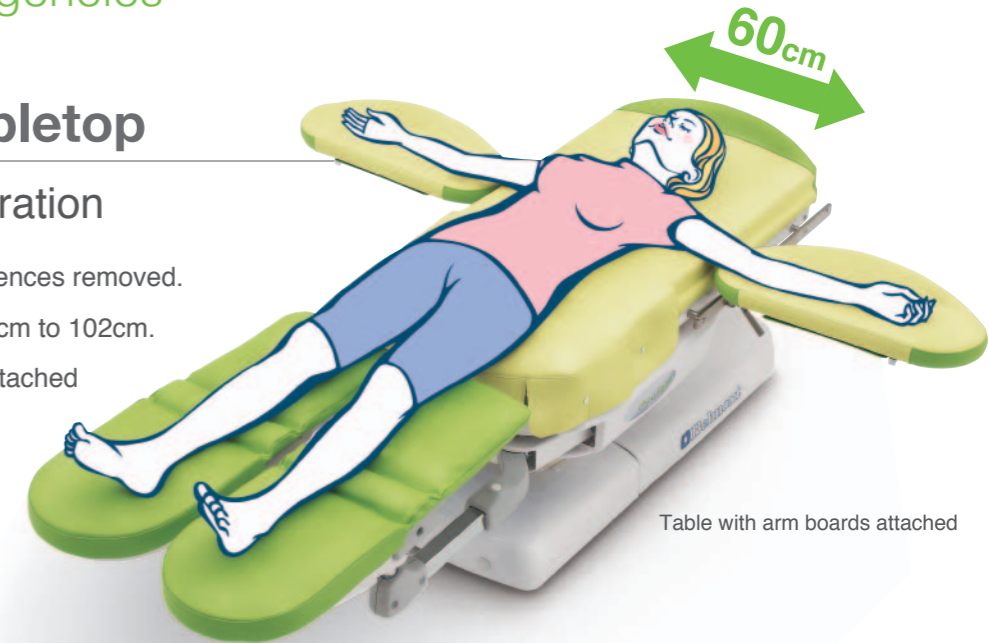


Table with arm boards attached

Trendelenburg and Lateral Tilting functions

Trendelenburg and lateral table top tilting by 15 deg. each

- Allows for securing a mother's position in case of emergencies including anesthetic procedures and responses to blood pressure reductions.
- Equipped with side rails for attaching necessary accessories.
- Offers a built-in battery.

Trendelenburg and Reverse Trendelenburg




Lateral Tilting



Smooth stretcher approaching

- Arm board can be detached to bring a stretcher alongside the table, facilitating mother transfer from/to the stretcher.

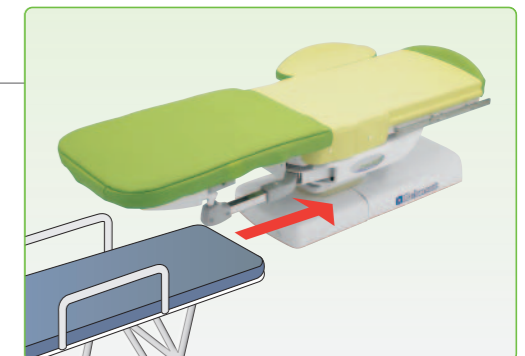

DELIVERIES

Various functions support mothers and practitioners during uneasy labor stage, delivery stage, and recovery stage in which mothers feel most at ease with their babies.


Improved comfort

Low resilient foam cushion of the mattress receives mother's body gently and firmly. This highly cushioned mattress on the table top reduces labor and delivery burdens on mothers and provides them with a sense of ease and comfort as if being enveloped in a warm hug.



Easy to get on/off the table

Initial height of 58cm enables mothers to easily get on and off the table in relief.

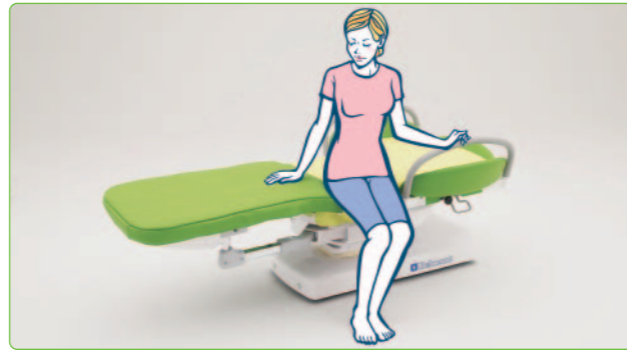

Provides mothers with a high comfort level and secure assistance, making mothers feel like being enveloped in a warm hug.



Electrically Operated Foot Support

Allows mothers to wide open their legs naturally with the ischial tuberosity as a point of support. Controllable with the foot switch system to easily adjust the opening angle of legs according to the delivery situation and operation.



Legs can be opened wide smoothly with ease.

Foot switch system enables practitioners to easily control the delivery table while assisting a mother.

The angle of the seat section can be independently adjusted according to the delivery situation and progress.

Various delivery positions supported



Supine Position

Sitting Lithotomy Position

Prone On All Fours Position

Lateral Position

Leg support for delivery in lateral position, which upholds mother's lower limb and supports mother's legs to help push in labor and delivery.

Hummocky leg support

Flexing action in a single operation allows the lower legs to be gently held on the heel rests.



Lower legs are held gently as if resting comfortably in a hammock.

Retractable newborn baby support board

Retractable type allows for approaching closer to the table during treatment and sutures. Slide mechanism supports delivery assistance from the lateral side.



When retracted



Support board is slid to secure space for assisting delivery.

Electro-hydraulic base lock mechanism


A combination of electro-hydraulic base lock mechanism and swivel casters allows the table to be easily moved, facilitating preparation and cleaning before and after delivery.



Knee Crutch:
DG-880-CN



Side rails:
DG-880-SR

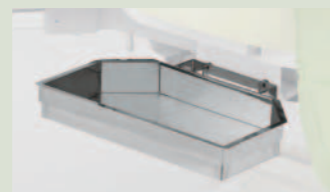

Douche Tray:
DG-880-XT (Picture)
Douche Tray with Drain:
DG-880-XR



Leg support for delivery in lateral position:
DG-880-SCG



Hand Switch:
DG-880-HSW



Strain Grip Above Head:
DG-880-HG



Cart for Accessories:
DG-880-OPC