

Obstetrics and Gynecology  
Examination Table

# DG-325

Obstetrics and Gynecology Treatment  
Tables for a Wide Range of Areas Including  
Consultation, Procedures, and Testing



Heel rest is an accessory sold separately.  
The product shown is DG-325 ES19 (Baby Pink).

## Color Selection



ES19 (Baby Pink)




ES9 (Ivory)



ES29 (Blue)

## Features

The Procedure Tray Is Equipped with a Sliding Function to Collect Body Fluids and Cleaning Solution that Splash Over a Wide Area.

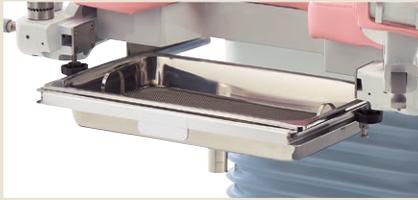

The Shape and Cushioning of the Upholstery Reduce the Burden of the Patient's Examination Position.



A Single Press of the Switch Moves the Table to the Programmed Setting or Return Position, Assisting efficient Medical Treatment.



The lowest positioning height makes it easy for elderly patient to get on and off the table.

The Electrohydraulic System Minimizes Vibration and Shock.



The Wireless Foot Control Eliminates the Stress Caused by Foot Control Wire in the Space Around the Table Base.



Infrared light transmitting through 360° allows the Foot Control to be positioned and angled in a wide range of locations during treatment.

An Auxiliary Switch is Provided on the Main Unit without the Use of a Foot Control.



P: Table Lifts to Programmed Height.  
R: Table Returns to Its Lowest Height.

## Option

Knee Crutches (one pair)\*<sup>1</sup>  
DG-325-C



Heel Supports (one pair) \*<sup>1</sup>  
DG-325-B



Arm Board\*<sup>2</sup>  
DG-325-A



Side Fence\*<sup>2</sup>  
DG-325-SFS



Sample image of the unit with the following optional items attached: knee crutches, side rails, and a side fence



\*<sup>1</sup>: Knee crutches and heel supports come with straps.  
\*<sup>2</sup>: Side rails (DG-325-A-SR) are required to install arm board and side fence. Side rails (DG-325-A-SR) are sold separately and are included for both sides in one unit.

Foot Board  
DG-325-H



Screen Frame  
UR-325-CBK




Instrument Holder  
(two pcs.) 1P06Z5A0




Cusco Holder  
1P06Z6A0



Sheet Rolls  
(12-roll pack)

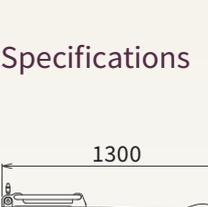


Drain Buckets  
Small (5 L), Large (10 L)

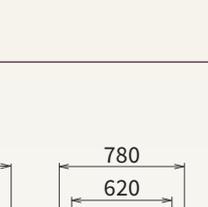



Plastic Seat Cover



## Specifications

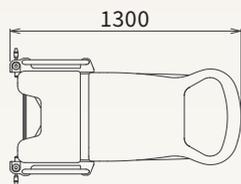

(Unit:mm)

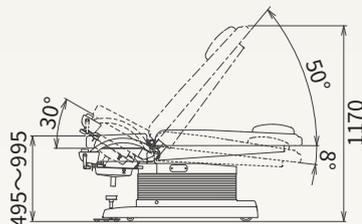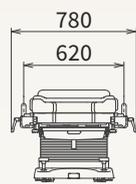

Seat Lifting Range	495 ~ 995 mm
Backrest Angle Range	-8° ~ 50°
Seat Angle Range	0° ~ 30°
Seat Lift to the Preset Position / Seat Return to the Lowest Position	Available / Available
Electrical Rating	AC 110 V, 50/60 Hz, 470/620 VA AC 120 V, 50/60 Hz, 510/680 VA AC 220 V, 50/60 Hz, 550/620 VA AC 230 V, 50/60 Hz, 510/580 VA AC 240 V, 50/60 Hz, 530/600 VA
Product Weight	135 kg

Operation mode

Non-continuous operation. Maximum operating time: 3 minutes, duty cycle: 1:5\*

\* Duty cycle is the ratio of downtime to continuous operating time of the hydraulic motor. (It varies with the ambient temperature.)

Standard accessories

• 1 Headrest Cushion • 1 Foot Control • 1 Drain Hose (with a band) • 1 Douche Tray (with wire netting) • 1 Plastic Seat Cover

**Belmont**

<https://medical.takarabelmont.co.jp/>

