

Urological Examination Table

UR-325

Urological Treatment Tables for a Wide Range of Areas Including Consultation, Procedures, and Testing



Features

Large Procedure Tray with 2 Stage Slide



The procedure tray is equipped with a mesh plate to prevent body fluids and cleaning solution that splash over a wide area.

The Shape and Cushioning of the Upholstery Reduce the Burden of the Patient's Examination Position.



A Single Press of the Switch Moves the Table to the Programmed Setting or Return Position, Assisting Efficient Medical Treatment.



The lowest positioning height makes it easy for elderly patient to get on and off the table.

The Electrohydraulic System Minimizes Vibration and Shock.



The Wireless Foot Control Eliminates the Stress Caused by Foot Control Wire in the Space Around the Table Base.



Infrared light transmitting through 360° allows the Foot Control to be positioned and angled in a wide range of locations during treatment.

An Auxiliary Switch is Provided on the Main Unit without the Use of a Foot Control.



P:
Table Lifts to
Programmed Height.

R:
Table Returns to
Its Lowest Height.

Option

Arm Board*
DG-325-A



Side Fence*
DG-325-SFS



Elbow Rests for the Doctor
UR-325-I



Foot Board
DG-325-H



Screen Frame
UR-325-CBK




Sheet Rolls
(12-roll pack)

Drain Buckets
Small (5 L), Large (10 L)

Plastic Seat Cover

* Side rails (DG-325-A-SR) are required to install arm board and side fence. Side rails (DG-325-A-SR) are sold separately and are included for both sides in one unit.

Specifications

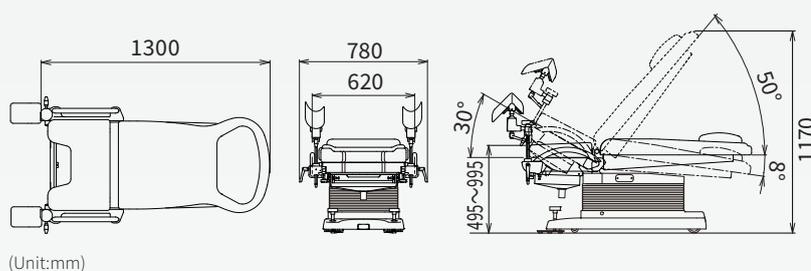

(Unit:mm)

Seat Lifting Range	495 ~ 995 mm
Backrest Angle Range	-8° ~ 50°
Seat Angle Range	0° ~ 30°
Seat Lift to the Preset Position / Seat Return to the Lowest Position	Available / Available
Electrical Rating	AC 110 V, 50/60 Hz, 470/620 VA AC 120 V, 50/60 Hz, 510/680 VA AC 220 V, 50/60 Hz, 550/620 VA AC 230 V, 50/60 Hz, 510/580 VA AC 240 V, 50/60 Hz, 530/600 VA
Product Weight	137 kg

Operation mode

Non-continuous operation. Maximum operating time: 3 minutes, duty cycle: 1:5*

* Duty cycle is the ratio of downtime to continuous operating time of the hydraulic motor. (It varies with the ambient temperature.)

Standard accessories

• 1 Headrest Cushion • 1 Foot Control • 1 Knee Rests (one pair) • 1 Douche Tray (with wire netting) • 1 Drain Hose (with a band) • 1 Plastic Seat Cover

Belmont

<https://medical.takarabelmont.co.jp/>

