


OPERATING TABLE
DR-2600 *Series*

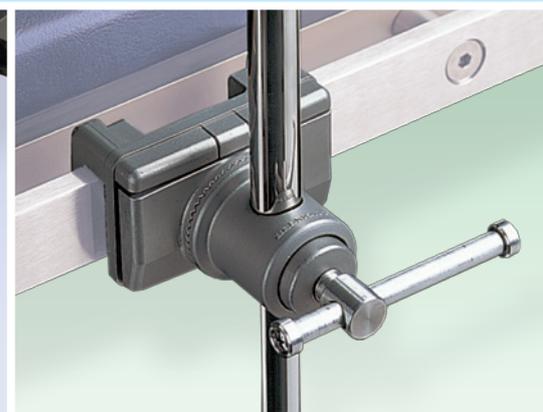

COMMON FEATURES OF DR-2600 SERIES

BELMONT'S RELIABLE ELECTRO-HYDRAULIC MECHANISM

High Power

Hydraulic device can generate greater power that takes little effort to use compared with electric motors.

Quiet and Smooth

The inertia of hydraulic device is so small and the movements keep quiet and smooth.

Long life

The lubricity and no rust of oil make the hydraulic device to require no parts replacement due to mechanical abrasion.

HORIZONTAL POSITION

LEVEL DETECTORS

Stops Reverse Trendelenburg, Trendelenburg, Lateral Tilt, and Back section (or Pelvic section) when horizontal during operation.

AUTO-RETURN-TO-LEVEL MECHANISM

A function to allow Reverse Trendelenburg, Trendelenburg, Lateral Tilt and Back section (or Pelvic section) to return to a horizontal position while pushing the button.

SWAPPING OF HEAD AND LEG SECTIONS

ONE table can be used for a variety of operations by swapping the positions of the head and leg sections.

Aiming for development hand in hand with world industry



When customers around the world are satisfied and industry prospers, we too develop and prosper. Based on this idea, we are constantly developing local business connections. Domestically and internationally, we work to provide the intangibles of satisfaction and comfort.




The pursuit of health is as old as Time. It continues—and will continue—to captivate people. Our wide range of products and Service are aimed at helping people like doctors and other medical staff succeed in this pursuit. We also help anyone who aspires to a healthy life. Belmont's products and services are in use when and where it counts most.



X-Ray FLUOROSCOPIC TABLE TOP

Full length cassette tunnel can hold the X-ray cassette. (sold separately)

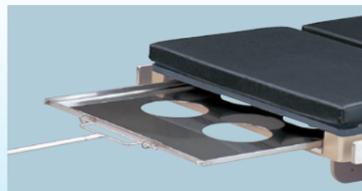

BASE COVER, MADE OF STAINLESS STEEL

High quality and well designed base

EMERGENCY STOP SWITCH



SAFETY CONSCIOUS 12VDC LOW-VOLTAGE HAND SWITCH CONTROLLER

The slim, lightweight remote control with icons ensures even greater operability, Erroneous operation is deterred with the "E" button, which needs to be depressed simultaneously with the desired function.



DR-2600
DR-2600-Y

DR-2600-S
DR-2600-Y-S

TABLE UNIT LOCK

Equipped with a base lock mechanism, controlled by our electro-hydraulic system, improving safety and operability.

Optional operating table accessories

Optional Belmont operating table accessories can easily be configured to appropriately support the position of the patient.



Arm Board



Anesthetic Screen Frame L-Shaped



Universal Lateral Support




Knee Crutches



Foot Control

DR-2600

Standard model

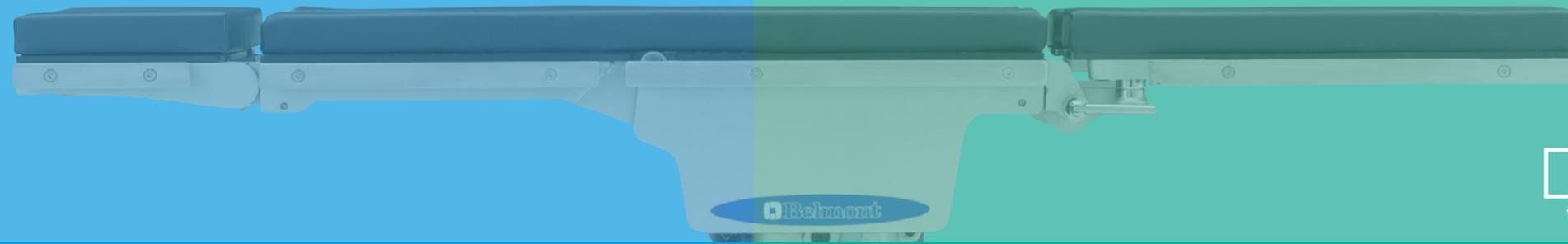

DR-2600-Y

With extra functionality

A cost-effective standard operating table designed for various surgical procedures.



- | | |
|--|---|
| <p>Table top</p> <ul style="list-style-type: none"> ● Swapping of Head and Leg section ● Full length cassette tunnel ● Consists of Head section, Back section, Pelvic section and Leg section <p>Other features</p> <ul style="list-style-type: none"> ● Level detectors ● Auto-Return-To-Level Mechanism | <p>Base</p> <ul style="list-style-type: none"> ● High quality stainless steel ● Emergency stop switch ● Table unit lock mechanism |
|--|---|

Batteries are mounted onto the existing standard model, with more convenient functions added to improve usability.



- | | |
|--|--|
| <p>Table top</p> <ul style="list-style-type: none"> ● Swapping of Head and Leg section ● Full length cassette tunnel ● Consists of Head section, Back section, Pelvic section and Leg section <p>Other features</p> <ul style="list-style-type: none"> ● Rechargeable battery ● Auto-Return-To-Level Mechanism ● Orientation Changing -Over Switch | <p>Base</p> <ul style="list-style-type: none"> ● High quality stainless steel ● Emergency stop switch ● Table unit lock mechanism ● Level detectors ● Override control panel |
|--|--|

BACK SECTION BENDING TYPE

Suitable for: Otorhinolaryngology, Ophthalmology, Dermatology, General surgery, etc.



Trendelenburg position



Kidney position

PELVIC SECTION BENDING TYPE [after swapping around Head section and Leg section]

Suitable for: Gynecology, Urology, General surgery, etc.



Supine position



Lithotomy position



Bent And Split Leg Plates position

Features

Main power, Battery Status Indicator

Orientation Changing-Over Switch

This function, after swapping the head and leg sections, ensures that the represented orientation for trendelenburg, Reverse Trendelenburg, and Lateral Tilt remain accurate on the remote control.



Override Control Panel

An override control panel is provided, allowing operation of all movements even when the remote control is not available.



DR-2600-S

Table top slide model



DR-2600-SY

Table top slide with extra functionality

Provision of a sliding function enables the area for fluoroscopy to be expanded to better cope with diverse operations.



- | | |
|---|---|
| <p>Table top</p> <ul style="list-style-type: none"> Swapping of Head and Leg section Full length cassette tunnel Consists of Head section, Back section, Pelvic section and Leg section 300mm table top slide <p>Other features</p> <ul style="list-style-type: none"> Level detectors Auto-Return-To-Level Mechanism | <p>Base</p> <ul style="list-style-type: none"> High quality stainless steel Emergency stop switch Table unit lock mechanism |
|---|---|

Table top longitudinal slide performed by electro-hydraulic mechanism



- | | |
|---|--|
| <p>Table top</p> <ul style="list-style-type: none"> Swapping of Head and Leg section Full length cassette tunnel Consists of Head section, Back section, Pelvic section and Leg section 300mm table top slide <p>Other features</p> <ul style="list-style-type: none"> Rechargeable battery Auto-Return-To-Level Mechanism Orientation Changing -Over Switch | <p>Base</p> <ul style="list-style-type: none"> High quality stainless steel Emergency stop switch Table unit lock mechanism Level detectors Override control panel |
|---|--|

COMMON FEATURES OF DR-2600-S & DR-2600-SY

POSITIONS



Kidney position



Trendelenburg position



Reverse Trendelenburg position




Lithotomy position




Sitting position


Back section bending mode (DR-2600-S / DR-2600-SY)

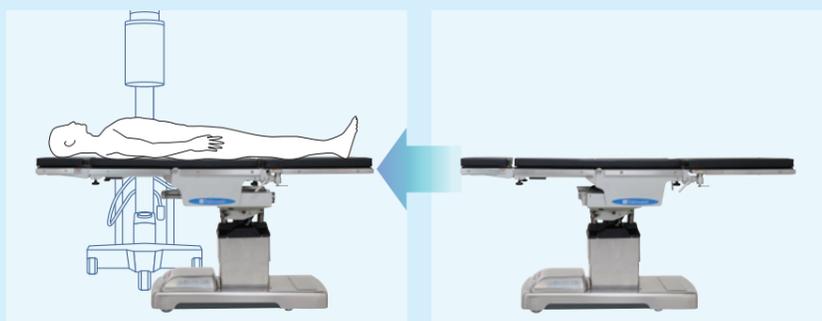

Fluoroscopic operation primarily for the torso

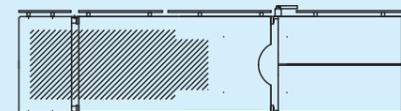

Pelvic section bending mode (DR-2600-S / DR-2600-SY)

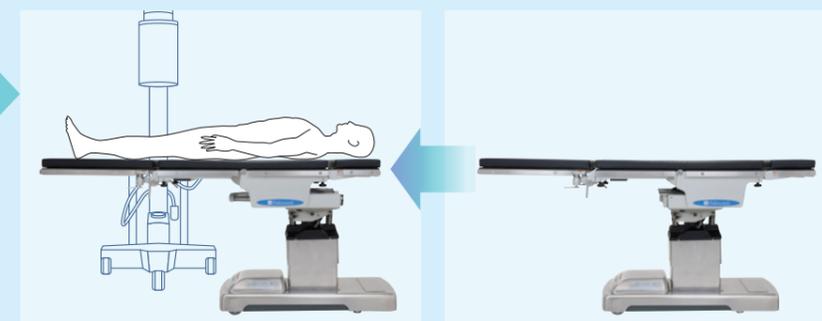

Fluoroscopic operation primarily for the lower limbs



Features

Main power, Battery Status Indicator


Orientation Changing-Over Switch

This function, after swapping the head and leg sections, ensures that the represented orientation for trendelenburg, Reverse Trendelenburg, and Lateral Tilt remain accurate on the remote control.

Override Control Panel

An override control panel is provided, allowing operation of all movements even when the remote control is not available.

TECHNICAL SPECIFICATIONS

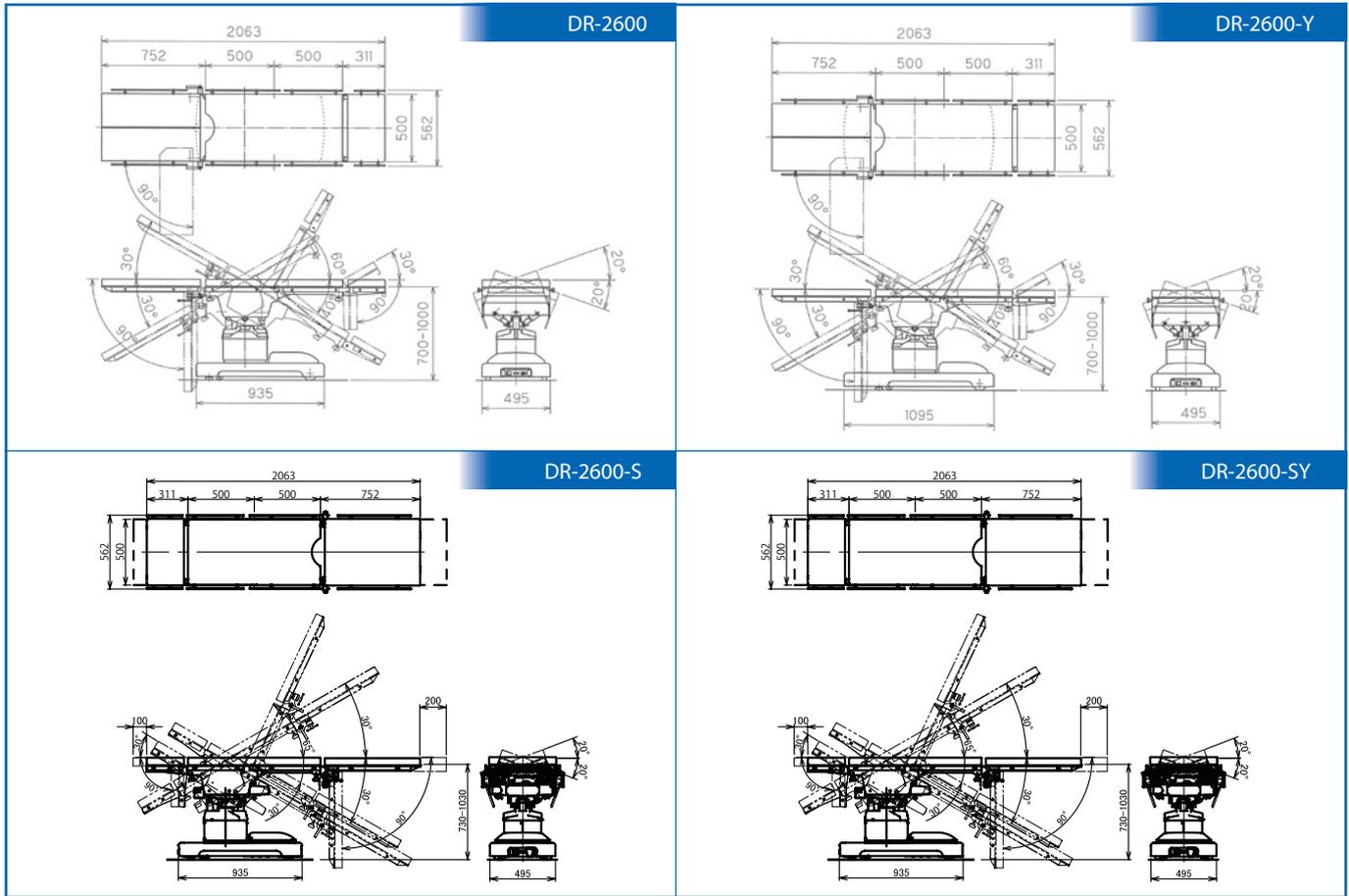

SPECIFICATIONS

		DR-2600	DR-2600-Y	DR-2600-S	DR-2600-SY
Table height	Electro-hydraulic mechanism	700-1000 mm	700-1000 mm	730-1030 mm	730-1030 mm
Trendenburg	30deg, Electro-hydraulic mechanism	○	○	○	○
Reverse trendelenburg	30deg, Electro-hydraulic mechanism	○	○	○	○
Lateral Tilt (right/ left)	20deg, Electro-hydraulic mechanism	○	○	○	○
Back section (or Pelvic)	up/ down, Electro-hydraulic mechanism	+60/-40 deg	+60/-40 deg	+65/-30 deg	+65/-30 deg
Table slide	head direction 200mm/ leg direction 100mm, Electro-hydraulic mechanism			○	○
Head section	up/ down, +30/-90 deg, Manual operation	○	○	○	○
Leg section bending	up/ down, 0/-90 deg, Manual operation	○	○	○	○
Leg section splitting	right/ left, 0/-90 deg, Manual operation	○	○	○	○
Auto-return-to-level mechanism	Trendelenburg/ Reverse trendelenburg, Lateraltilt, and Back section(or Pelvic)	○	○	○	○
Level detectors	Trendelenburg/ Reverse trendelenburg, Lateraltilt, and Back section(or Pelvic)	○	○	○	○
Override control panel			○		○
Rechargeable Battery			○		○
Orientation changing over switch					○
Rated voltage	AC. 110/ 120/ 220/ 230/ 240V	○	○	○	○
Operating voltage	DC-12V	○	○	○	○
Product weight	DR-2600 220KG/ DR-2600-Y 240KG/ DR-2600-S 250KG/ DR-2600SY 270KG				

The manufacturer reserves the right to alter specifications without prior notice.



TAKARA BELMONT CORPORATION
 2-1-1, Higashishinsaibashi,
 Chuo-ku, Osaka, 542-0083, Japan
 TEL : 81-6-6213-5945
 FAX : 81-6-6212-3680