

CR-LPR

Series

MOBILE EXAMINATION
&
TREATMENT TABLE




Concept

Examination, Treatment, and Patient Transport

—all with just this one table

A One-Stop Solution That Improves Workflow Efficiency Throughout the Clinical Environment

We propose a mobile table that combines full examination & treatment table functionality with excellent mobility for patient transport to other areas, including recovery rooms.

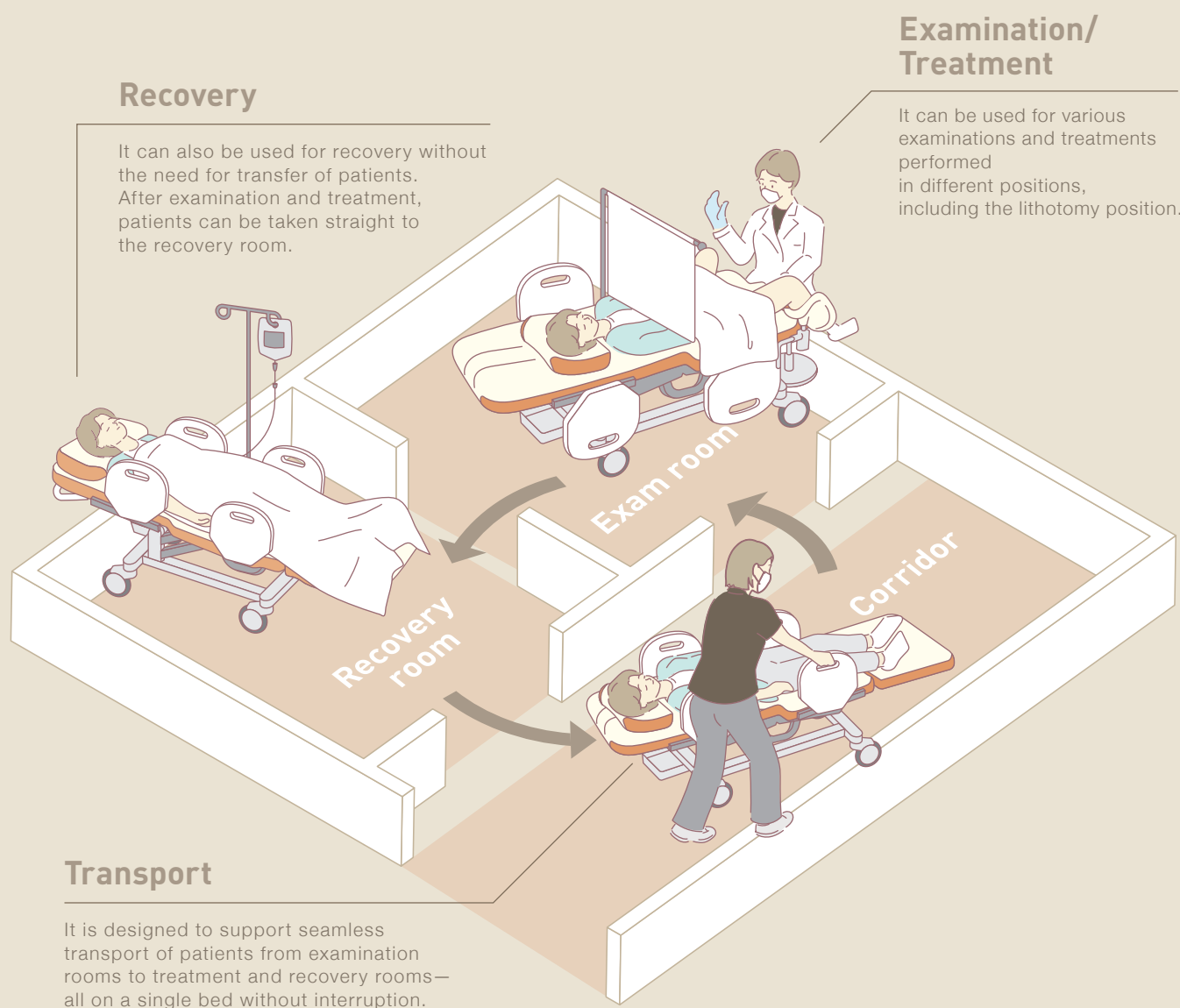

M type Multi Type [CR-LPRM]

Adaptable to a Wide Range of Clinical Scenarios

Thanks to its compatibility with variety of accessories, it supports a wide range of patient positioning, leveraging BELMONT's expertise as a surgical table manufacturer.




It enables a variety of patient positions that accommodate various clinical settings. By selecting the most appropriate position for each situation, you can use it widely across various medical specialties.

* The product shown is M type.
The knee crutches(CR-LPR-CZ) is not compatible with C type.

Quick Reclining Operation


Effortless operation supported by a gas cylinder regardless of a patient's physical build. The slip-resistant mechanism that lowers the seat section as the backrest is raised prevents the patient from sliding.

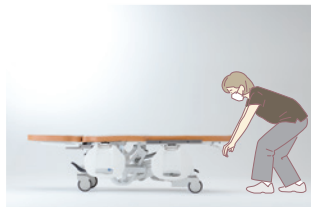

C type Core Type [CR-LPRC]

Maximizing Space Efficiency

The linked leg-rest bending mechanism works in conjunction with the backrest movement, allowing the chair to save space when stored in the standby position.

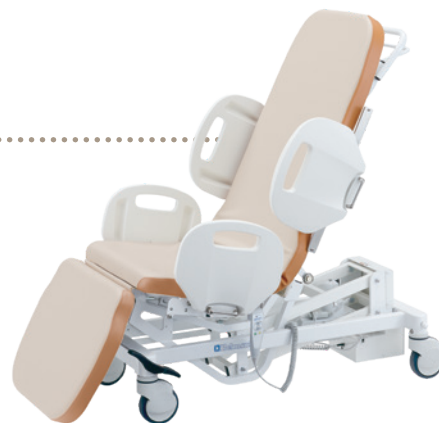

Features



Better Mobility with Durable Castors

Single-pedal operation enables switching between four-wheel free, directional-travel, and four-wheel lock. When totally locked, castors do not rotate or swivel, ensuring stability during examination and treatments.



Lithium-ion Battery

It uses a compact lithium-ion battery with a long runtime and durability for repeated recharging compared to lead-acid batteries.

The table can be moved up and down approximately 100 times on a single charge (when new). The audible alarm and LED indicator will alert you when the battery level is low.



Quick and Effortless Setup


It is designed to support emergency transport and side-by-side stretcher positioning.



Rotatable Side Fence

The rotatable side fences do not require space when switching between the use and storage positions.

When not in use, they can be stored below the tabletop, which provides easy accessibility during treatments or when transferring the patient onto or off of the table.

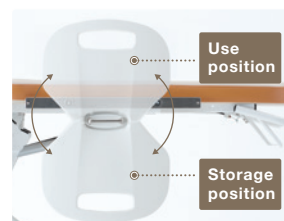

One-Handed Operation

Each fence can be operated with one hand. When moved to the use or storage position, it is automatically locked by releasing it. (It cannot be unlocked even when pushed from inside.)





Lithotomy Position

It accommodates the lithotomy position, making it suitable for examinations and treatments in gynecology and urology applications.



* The product shown is M type.
The knee crutches(CR-LPR-CZ) is not compatible with C type.

Height Adjustable Anytime, Anywhere

Height adjustment available throughout the transport process and at the destination. Since it is battery-powered, no power cord is required.



Hand Controller with Two Memory Settings

With two preset switches, you can adjust the table to your preferred height with a single press.

OPTION



Foot Controller with Two Memory Settings

The table can also be operated with an optional foot controller.



Stand CPR

High-rigidity construction engineered for demanding emergency scenarios.

Accessories Customization

Common Accessories



1
EB-BFT-M
HEADREST CUSHION (M)



2
CR-LPR-M
HEADREST CUSHION (L)



3
EX-CS-AT
ARM BOARD



4
CR-LPR-ASR
ARM SUPPORT R




5
DR-OPTM-G05
RAISED ARM BOARD



6
DR-OPTM-SIBK/SAIK
TRANSLUCENT ARMBOARD
* Also available without a cushion.
(DR-OPTM-SIA/SIB)



7
EX-CS-E
LATERAL SUPPORT

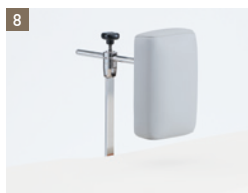


8
EX-CS-EN
UNIVERSAL LATERAL SUPPORT



9
CR-LPR-BH
OXYGEN CYLINDER HOLDER
* Holder Dimensions
-Length : 440 mm -Inside Diameter : 109 mm
* Oxygen cylinder is not included.



10
EX-CS-BLM
L-SHAPED SCREEN FRAME

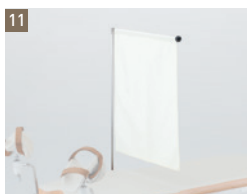

11
CR-LPR-CB
SCREEN FRAME



12
EX-CS-TD
L-SHAPED TRAY HOLDER



13
CR-LPR-H
LOWER LIMB SUPPORT
* Side Rail (CR-LPRC-SR, 1 set, one side)
required for mounting the CR-LPR-H on
the CR-LPRC.



14
CR-LPR-KM
TRAY



15
CR-LPR-FK
FOOT CTRL WITH TRAY



16
EB-BFT-SC
SIDE FENCE COVER



17
CR-LPR-B
BODY RESTRAINT STRAP



18
CR-LPR-DC/CR-LPR-DC-BK
IV POLE FOR CENTER*
/ BRACKET ONLY
* The bracket for 1pc included.



19
DR-OPTM-ZNN
IV POLE



20
CR-LPR-PH/DG-OPT-R35
ROLL SHEETS HOLDER/
DISPOSABLE ROLL SHEETS



21
CR-LPR-AC
ACCESSORIES CART



22
DR-LIBP
LITHIUM-ION BATTERY
* Battery charger is not included.




23
DR-LICHG
BATTERY CHARGER
FOR LITHIUM-ION BATTERY



24
CASTER LOCK PEDAL

Lithotomy Position [M type]

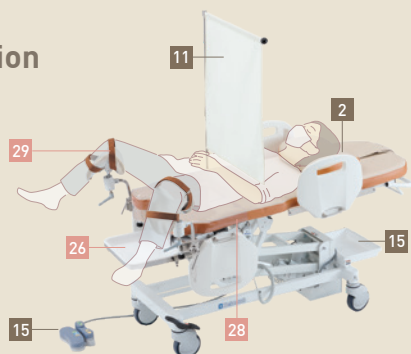


Supine Position [M type]




*The Screen is not included.


Lateral Position [C type]

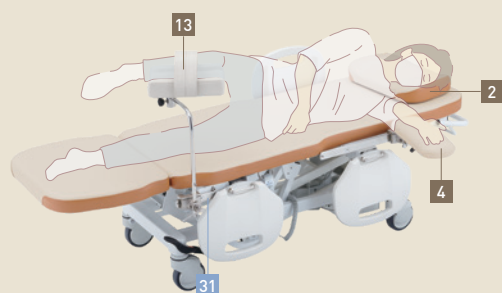

Seat Position [C type]

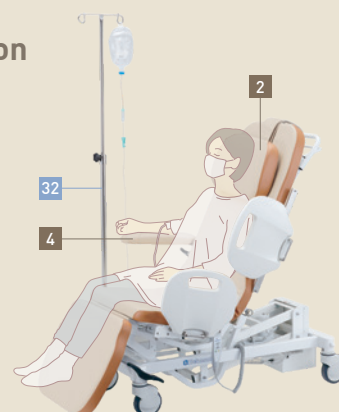

For M-Type Only

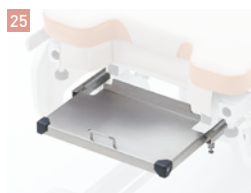

CR-LPR-BT
INSTRUMENT TRAY HOLDER



CR-LPR-XT/CR-LPR-XTR
PROCEDURE TRAY/
PROCEDURE TRAY (with Drain Hose)



CUSCO HOLDER
Procedure Tray (CR-LPR-XT/CR-LPR-XTR,
1 set) required for mounting the Cusco
Holder.



CR-LPR-SC
SEAT COVER



CR-LPR-CZ
KNEE CRUTCHES (S)



CR-LPR-CN
KNEE CRUTCHES (L)


For C-Type Only



CR-LPRC-SR
SIDERAIL FOR LPRC
* 1Set(left and right)



CR-LPRC-DC/CR-LPRC-DC-BK
IV POLE FOR FOOT DIRECTION*
/BRACKET ONLY
* The bracket for 1pc included.

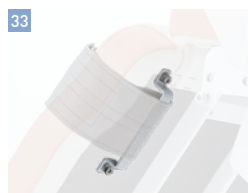

CR-LPRC-BL
HOLDER FOR BODY RESTRAINT
STRAP

* Body restraint strap is not included.
1Set (left and right)

*Options not included in this catalogue cannot be installed.

Color Selection



ES9(Ivory)
×ES14(Camel)



ES9(Ivory)
×ES20(Coral Pink)



ES09(Ivory)
×ES23(Aldo Blue)

* The above colors are available for both M/C types.

Antibacterial performance and Alcohol resistance



Antibacterial performance

Tested in accordance with JIS Z 2801:2012 Antimicrobial products - Test for antimicrobial activity and efficacy (Test bacteria: Staphylococcus aureus and Escherichia coli)



Alcohol resistance

Tested in accordance with our standard protocol


Specifications

Model	CR-LPRM (M type), CR-LPRC (C type)
Dimensions	600mm×1850mm (tabletop)
Table height range	500mm~930mm (electric motor)
Side fence	Rotatable & retractable
Castor	Four-wheel free, Directional-travel, Total lock


Electrical rating	DC21.6V 155VA
Product weight	93kg
Max. load capacity	180kg
Standard accessories	Hand controller, Battery pack, Charger

Dimensions

[M type]

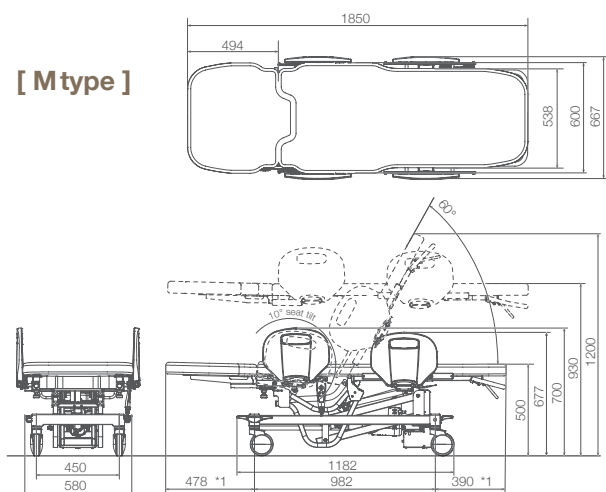

[C type]

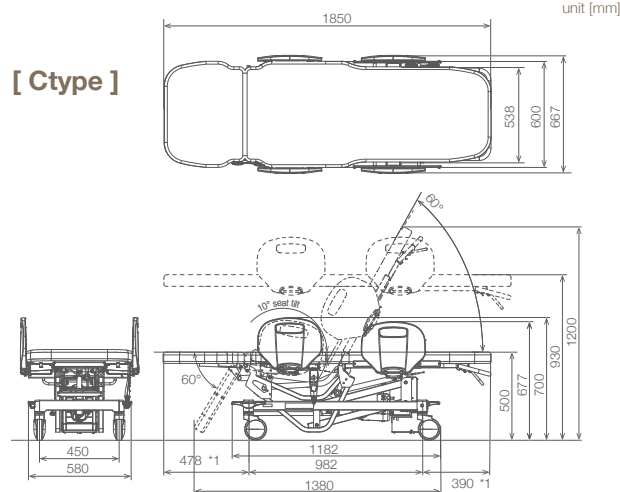

*1 Measured at the lowest position. Dimensions vary depending on the device height.



<https://megb.takarabelmont.co.jp/>



TAKARA BELMONT CORPORATION 2-1-1, HIGASHISHINSAIBASHI, CHUO-KU, OSAKA, 542-0083, JAPAN +81-6-6213-5945

The specifications and colors are subject to change without notice.

REF M-034 2607A MI